


## Aspect100RT Multi-axis

### PRODUCT DATA

#### Features:

- Low attenuation, carbon fiber table construction
- Tabletop: 24" x 80"
- Imaging area: 24" x 63"
- Lateral roll:  $\pm 20^\circ$
- Trendelenburg:  $\pm 15^\circ$
- Height range: 30"- 40"
- Capacity: 450#
- Table weight: ~325#
- Motorized elevation
- Locking swivel casters
- Low attenuation table pad
- Foot switch and hand control with return to level function

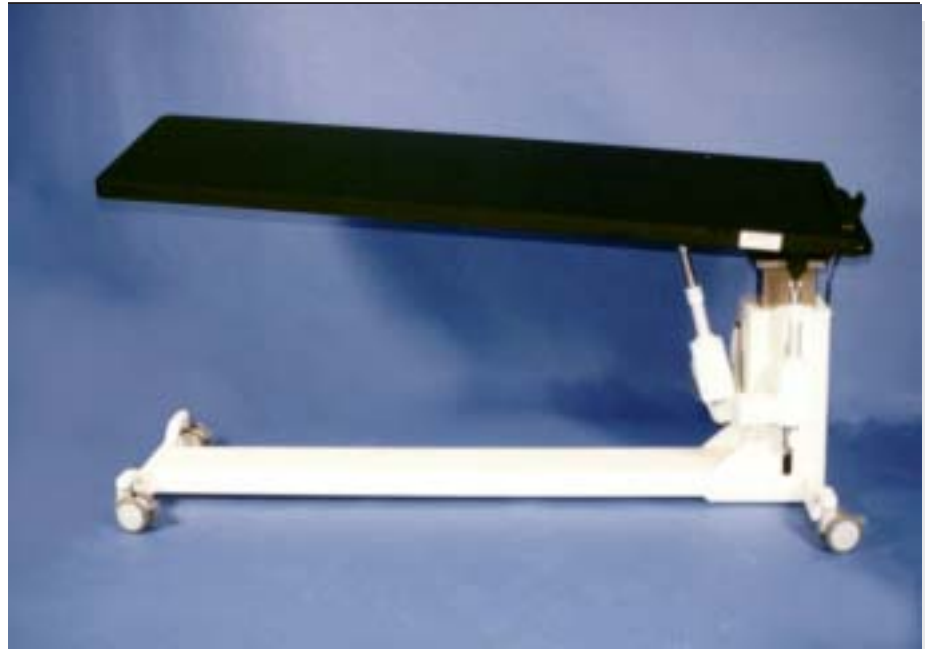

### Aspect100RT Fluoroscopic Examining Table

Image Diagnostics' *Aspect100RT* is specifically designed for ease of use and high quality image outcome. *Aspect* tables benefit from a trim, lightweight design, ensuring a single operator can position the table with ease. Features unique to the *Aspect* table series provide the ultimate in table maneuverability, operator access and patient comfort. The Image Diagnostics *Aspect100RT* is designed with full surgical articulation for use in vascular, endoscopic, pain management and orthopedic procedures requiring multi-axis positioning.

The innovative, cantilevered table design of the *Aspect100RT* significantly enhances mobile C-arm positioning and provides visualization capabilities free from obstructions or artifacts. A unique, center-rib base design allows the operator to stand tableside without having to step on or around casters. Approaches to the table can be made from either side as well as the head-end, allowing for enhanced patient access with little or no need for repositioning.

The *Aspect100RT* is designed with an extended base that eliminates the need for heavy ballast for stability and helps protect the tabletop from impacts. Vertical travel capabilities of the *Aspect100RT* allow patient placement at a comfortable working height and facilitate its successful integration into a variety of mobile C-arm procedures.

A variety of table accessories are available for patient positioning and operator radiation protection.

# Aspect100RT Fluoroscopic Examining Table

---

## FEATURES AND BENEFITS

- Carbon fiber tabletop is 80" in length with 62.5" imaging area; for most patients, provides total body coverage with minimal need for C-arm repositioning.
- Tabletop is constructed from low attenuation carbon fiber; enhances image quality while reducing dose.
- Table width is 24"; provides support for larger patients while maintaining excellent access for lateral C-arm positioning.
- Table articulation: Trendelenburg  $\pm 15^\circ$  and lateral roll of  $\pm 20^\circ$  allows optimal patient positioning for working ease and quality imaging.
- Motorized elevation with height range of 30" to 40" provides one of the lowest available retracted heights; facilitates patient loading and unloading.
- One button "return to level" allows operator to return to a home position and reposition table top for loading and unloading.

*For a complete listing of Image Diagnostics' tables and products, visit our website at [www.imagediagnostics.com](http://www.imagediagnostics.com).*



*All Aspect series tables and accessories are backed by a three year warranty on parts.  
Image Diagnostics, Inc. reserves the right to update product specifications, without notification, as required.  
[www.imagediagnostics.com](http://www.imagediagnostics.com)*

100RT